

This product is **for research use only** (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

## Product no AS10 1201

Donkey anti-Mouse IgG (H&L), DyLight® 488 conjugated, min, cross-reactivity to bovine, chicken, goat, guinea pig, hamster, horse, human, rabbit, rat or sheep IgG

## **Product information**

**Immunogen** Purified Mouse IgG, whole molecule

**Host** Donkey

Clonality Polyclonal

Purity Immunogen affinity purified donkey IgG.

Format Lyophilized

Quantity 1 mg

Reconstitution

For reconstitution add 1.1 ml of sterile water. Let it stand 30 minutes at room temperature to dissolve. Prepare fresh working dilutions daily

Storage

Store lyophilized material at 2-8 °C. Product is stable for 4 weeks at 2-8 °C after rehydration. For long time storage after reconstitution, dilute the antibody solution with glycerol to a final concentration of 50% glycerol and store as liquid at -20 °C, to prevent loss of enzymatic activity. For example, if you have reconstituted 1 mg of antibody in 1.1 ml of sterile water add 1.1 ml of glycerol. Such solution will not freeze in -20°C, If you are using a 1:5000 dilution prior to diluting with glycerol, then you would need to use a 1:2500 dilution after adding glycerol. Prepare working dilution prior to use and then discard. Be sure to mix well but without foaming.

**Additional information** 

Conjugate is present in 10 mM Sodium Phosphate, 0.15 M Sodium Chloride, pH 7.2, 1 % (w/v) BSA, Protease/lgG free. 0.05 % (w/v) sodium azide is added as preservative.

Based on immunoelectrophoresis, this antibody reacts with: heavy ( ) chains on mouse IgG, light chains on all mouse immunoglobulins


No reactivity is observed to: non-immunoglobulin mouse serum proteins, IgG from bovine, chicken, goat, guinea pig, hamster, horse, human, rabbit, rat or sheep

## Application information

**Recommended dilution** 1:20-1:2000 for most applications

Selected references

Maurya et al. (2017). Hedgehog signaling regulates ciliary localization of mouse odorant receptors. Proc Natl Acad Sci USA. 2017 Oct 31;114(44):E9386-E9394. doi: 10.1073/pnas.1708321114.

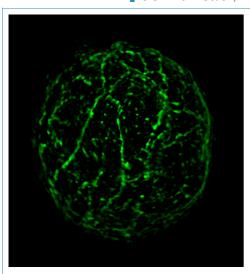

Sample: Maize protoplasts fixed in 4% paraformaldehyde in 1X PBS pH 7.4 for 1 hour and washed 3 times in 1XPBS buffer

Primary antibodies: Agrisera AS21 4615 actin monoclonal; clone 14H4G8;

Secondary antibodies: Donkey anti-Mouse IgG (H&L), DyLight® 488 conjugated AS10 1201 (Agrisera)



This product is for research use only (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

Other reagents: 30% Bovine Serum Albumin solution, Sigma Aldrich Cat # A8577-50ML; Phosphate Buffered Saline 10X, Electron Microscopy Sciences Cat # 19342-10; Imaging spacer, Millipore Sigma Grace Bio-Labs SecureSeal™ imaging spacer Cat # GBL654008

## Protocol:

- 1. Block isolated protoplasts in a blocking buffer (2% BSA in 1x PBS) for 1 h at room temperature.
- 2. Dilute the primary antibody to 1:500 in a blocking buffer, mix well, and spin down at 150 x g for 3 min:
- 3. Spin down protoplasts in a blocking buffer at 150 g for 1 -3 min and remove as much supernatant as possible using a pipette tip with a cut end (wide opening).
- 4. Resuspend protoplasts in diluted antibody and incubate for 1 h at room temperature.
- 5. Spin down protoplasts at 150 g for 1 min and remove as much supernatant as possible using a pipette tip with a cut end.
- 6. Wash protoplasts in a blocking buffer 3 times for 5 minutes each. Spin and resuspend as described in step 5 between each washing step.
- 7. Dilute the secondary antibody AS10 1201 to 1: 400, mix well, and spin down at 150 x g for 3 min:
- 8. Incubate protoplasts in diluted antibody for 1 h at room temperature.
- 9. Repeat the washing steps as described in 5-6.
- 10. Resuspend a protoplast pellet in a small volume of PBS 1X, pH 7.4.
- 11. Prepare microscopy slide: attach imaging spacer on top of the slide, peal off top adhesive membrane, add protoplast suspension on to the slide inside the spacer opening, cover with 22x22 mm cover glass No 1.5.
- 12. Imaging: Zeiss Elyra 7 SIM

Courtesy of Dr. Kirk Czymmek, The Danforth Plant Science Centre, USA