

This product is for research use only (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

Product no AS16 3155

## Anti-ADK | Adenylate kinase

## **Product information**

Immunogen Protein derived from Escherichia coli, Uniprot: P69441

**Host** Rabbit

Clonality Polyclonal

**Purity** Serum

Format Lyophilized

Quantity 50 ul

**Reconstitution** For reconstitution add 50 μl, of sterile water

Storage | Store lyophilized/reconstituted at -20°C; once reconstituted make aliquots to avoid repeated freeze-thaw cycles. Please remember to spin the tubes briefly prior to opening them to avoid any losses that might occur from material adhering to

the cap or sides of the tube.

## Application information

Recommended dilution 1: 100 000 (WB)

Expected | apparent

MW

23.5 kDa (E.coli), 24.3 kDa (S.cerevisiae)

Confirmed reactivity Escherichia coli

Predicted reactivity


Not reactive in No confirmed exceptions from predicted reactivity are currently known

Selected references

Tükenmez et al. (2016). Linkage between Fitness of Yeast Cells and Adenylate Kinase Catalysis. PLoS One. 2016 Sep

19;11(9):e0163115. doi: 10.1371/journal.pone.0163115. eCollection 2016.

## Application example

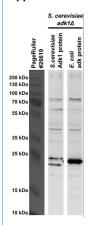

Protein samples were prepared using TCA-extraction method. S. cerevisiae adk1Δ strains expressing either S. cerevisiae Adk1 protein (24.3 kDa) or E. coli adk protein (23.5 kDa) were cultivated logarithmically, and 5 OD-units of cells were collected for TCA-extraction method. In the end. protein pellets were re-suspended with 300 μl of TCA-Laemmli Loading Buffer and denatured at 100°C for 10 minutes. Total protein samples (8 μl of each sample) were separated by 12% SDS-PAGE and transferred to nitrocellulose membranes using a wet transfer protocol. The membranes were blocked with 5% dry-milk in 1x TBST for 1.5 hours at room temperature (RT) with agitation and then they were incubated with primary antibody (anti-adk antibody, 1:2,000) in 5% dry-milk in 1x TBST overnight at 4°C with agitation. The primary antibody solution was decanted and the membranes were washed 3 times for 10 minutes with 1x TBST with agitation. The membranes were then incubated with secondary antibody (HRP conjugated anti-rabbit IgG, 1:4,000) in 5% dry-milk in 1x TBST for 1 hour at RT with agitation. The membranes were washed as above and the signals were detected using chemiluminescent detection reagent and exposing the membranes in LAS4000 machine (GE-healthcare) for 10 seconds. Followed-up experiments revealed that lowering the primary antibody concentration from 1:2,000 to 1:50,000-100,000 reduces the unspecific background signals and lowers the risk of target signal oversaturation.

Courtesy of Dr. Hasan Tukenmez, Umeå University, Sweden