

This product is for research use only (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

Product no AS16 3935

Anti-E(z) | Histone-lysine N-methyltransferase E(z)

Product information

Immunogen GST-fusion with amino acids 8-190 of the *Drosophila melanogaster* E(z) protein, UniProt: P42124

Host Rabbit

Clonality Polyclonal

Purity Immunogen affinity purified and depleted of anti-GST antibodies, in PBS pH 7.4.

Format Lyophilized

Quantity 50 μg

Reconstitution For reconstitution add 50 μl of sterile water

Storage Lyophilized antibody can be stored at -20 °C for up to 3 years. Re-constituted antibody can be stored at 4 °C for several days to weeks. Once reconstituted make aliquots to avoid repeated freeze-thaw cycles. Please remember to spin the tubes briefly prior to opening them to avoid any losses that might occur from material adhering to the cap or sides of the

Application information

Recommended dilution 3 μg/IP (ChIP)

Expected | apparent MW 86.9 kDa

Confirmed reactivity | Drosophila melanogaster

Predicted reactivity Bactrocera cucurbitae, Bactrocera dorsalis, Bactrocera latifrons, Ceratitis capitata, Lucilia cuprina

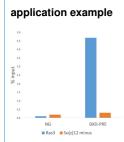


Figure 1. ChIP recovery.

ChIP and qPCR analysis were done as described [Schwartz YB, Kahn TG, Nix DA, Li XY, Bourgon R, et al. (2006) Genome-wide analysis of Polycomb targets in Drosophila melanogaster. Nat Genet 38: 700–705. doi: 10.1038/ng1817]. Chromatin from Su(z)12 mutant cell line (Kahn T. et all, submitted) served as a negative control, chromatin from Ras3 cells served as a positive control. Quantitative PCR was performed with primers specific for BXD-PRE of Ubx gene (Polycomb target gene in repressed state), used as positive controls, and for intergenic region, used as negative control. Figure 1 shows the ChIP recovery, measured by qPCR as a % of input (the relative amount of immunoprecipitated DNA compared to input DNA). Chromatin from 5x107 cells and 3mg of anti-E(z) antibody were used for each ChIP reaction.

Courtesy of Dr. Tatyana Khan, Umeå University, Sweden.