

This product is for research use only (not for diagnostic or therapeutic use)

contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

## Product no AS23 4895

### Anti-PTOX | Plastid terminal oxidase

#### **Product information**

Immunogen Recombinant PTOX of Arabidopsis thaliana, expressed in E.coli, UniProt: Q56X5, TAIR: AT4G22260

**Host** Rabbit

Clonality Polyclonal

**Purity** Antigen affinity purified serum, in PBS pH 7.4

Format Lyophilized

Quantity 50 μg

**Reconstitution** For reconstitution add 50 μl, of sterile or deionized water.

Store lyophilized/reconstituted at -20°C; once reconstituted make aliquots to avoid repeated freeze-thaw cycles. Please, remember to spin tubes briefly prior to opening them to avoid any losses that might occur from lyophilized

material adhering to the cap or sides of the tubes.

# Application information

**Recommended dilution** 1 : 2000 - 1: 5000 (WB)

Expected | apparent

40.5 | 34.3 kDa (due to N-terminal or C-terminal processing)


**Confirmed reactivity** Arabidopsis thaliana, Chlamydomonas reinhardtii

Predicted reactivity

Camelina sativa, Brassica napus, Capsella rubella, Capsicum annuum, Citrus sinensis, Cucumis sativus, Eutrema salsugineumm, Glycine max, Gossypium hirsutum, Malus domestica, Medicago truncalula, Nicotiana tabacum, Oryza sativa Japonica Group, Physcomitrium patens, Populus alba x Populus x berolinensis, Pisum sativum, Pistacia vera, Ricinus communis, Solanum lycopersicum, Spinacia oleracea, Sorghum bicolor, Vitis vinifera

Species of your interest not listed? Contact us

**Selected references** To be added when available, antibody available in May 2024.

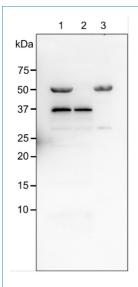

- 1 Chloroplast samples corresponding to 3.0 µg chlorophyll of Arabidopsis thaliana
- 2 Thylakoid membrane samples corresponding to 3.0 μg chlorophyll of Arabidopsis thaliana
- 3 Stromal fraction samples corresponding to 3.0 µg chlorophyll of Arabidopsis thaliana

Left panel: MW markers

Protein samples corresponding to 1.0 µg chlorophyll of chloroplasts extracted from Arabidopsis thaliana. Chloroplast samples were isolated with 20 mM Tricine-NaOH, pH 8.4, containing 400 mM sorbitol, 5 mM MgCl<sub>2</sub>, 5 mM MnCl<sub>2</sub>, 2 mM EDTA, 10 mM NaHCO<sub>3</sub>, 0.5% (w/v) bovine serum albumin (BSA) and 5 mM ascorbate: and denatured with exact buffer (62.5 mM Tris-HCl, pH 6.8, containing 2% (w/v) SDS, 10% (w/v) glycerol, 0.005% (w/v) bromophenol blue, and 50 mM dithiothreitol) at 95 °C/5 min. Samples were separated at RT on 15% acrylamide gel and blotted for



This product is for research use only (not for diagnostic or therapeutic use)


contact: support@agrisera.com

Agrisera AB | Box 57 | SE-91121 Vännäs | Sweden | +46 (0)935 33 000 | www.agrisera.com

0.5 h to PVDF (pore size of 0.2  $\mu$ m), using: semi-dry transfer at RT. Blot was blocked with 5 % milk for: 1h/RT with agitation. The blot was incubated in the primary antibody at a dilution of 1: 2,000 for 1h/RT with agitation in TBS-T. The antibody solution was decanted, and the blot was rinsed briefly twice, then washed 3 times for 5 min in TBS-T at RT with agitation. Blot was incubated in matching secondary antibody (<u>AS09 602</u> anti-rabbit IgG horse radish peroxidase conjugated) diluted to 1: 25,000 in for 1h/RT with agitation. The blot was washed as above and developed with a following chemiluminescent detection reagent: <u>AS16 ECL-N-10</u> AgriseraBright (mid picogram). Exposure time was 3 minutes.

Note: The band at 50 kDa is a cross reactivity with Rubisco.

Courtesy of Prof. Yuki Okegawa, Okayama University, Japan

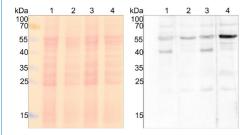

#### Samples:

- 11 Chlamydomonas reinhardtii, wild type cc124 (mt -)
- 2 Chlamydomonas reinhardtii ptox2.34 (mt -)
- 3 Chlamydomonas reinhardtii wild type cc125 (mt +)
- 4 Chlamydomonas reinhardtii ptox2.24 (mt +)

1x10e6 cells per lane of total protein extracted with Chloroform/Methanol protein precipitation protocol (dx.doi.org/10.17504/protocols.io.icecate) from *Chlamydomonas reinhardtii*. Exact buffer components were: Count *Chlamydomonas reinhardtii* cell density in TAP cultures (early log phase), spin down 15 ml culture (4000 g, 5 min, 22 °C) and resuspend pellet in 1 ml supernatant. Use 200 μl resuspended cells for protein precipitation in 2 ml Eppendorf tubes at room temperature (dx.doi.org/10.17504/protocols.io.icecate). Pellets were dried from MeOH and resuspended in 1x SDS PAGE running buffer (50 mM Tris/HCl pH 6.8, 50 mM DTT, 50 mM Na<sub>2</sub>CO<sub>3</sub>, 2% (w/v) SDS, 12% (v/v) Glycerol, 0.01% (w/v) Bromphenolblau, 0.01% (w/v) Servablau G) so that 20 μl corresponds to 1x10e6 cells. Sample was denatured (20 min, 22 °C) and spun down before loading. 13% SDS-PAGE Laemmli 1 h to nitrocellulose (pore size of 0.45 μm) semi-dry at room temperature. Ponceau S staining of total proteins to check loading and transfer efficiency. Follow recommended protocol for the remaining procedure (use agitation throughout; no double rinsing but in between blocking & antibody changes include 3x washing steps with TBS-T, 5 min/RT). Blocking: 5% non-fat milk 1h/RT, primary anti-PTOX antibody (Product No AS23 4895, Lot 2309) 1:1000 in TBS-T ON/4C, secondary antibody (AS09 602-trial) 1:25000 1h/RT, detection first with Agrisera Bright <u>AS16 ECL-N</u>. Exposure time was 3 s.

Courtesy of Dr. Felix Buchert, University of Munster, Germany